

Easy DRY™ Gel Drying System

Cost-Effective and Crack-Free Reversible Gel Drying


Product Highlights

- Dry gels for permanent records
- No gel drying equipment required
- Proprietary gel drying solution prevents cracking of gels
- Comes with pre-cut gel drying sheets (50)
- Comes with frame and rust-free clamps
- Most economic and user-friendly system for drying gel



Comparison

Proteins were loaded on 12% SDS-TG gels. Following electrophoresis, proteins were stained with NRL Ultra-FAST Coomassie stain or Ultra-FAST Silver Stain. Stained gels were dried with NRL Easy-GEL drying system with or without Gel Drying Solution (GDS-004). Gels dried with GDS-004 had no cracking

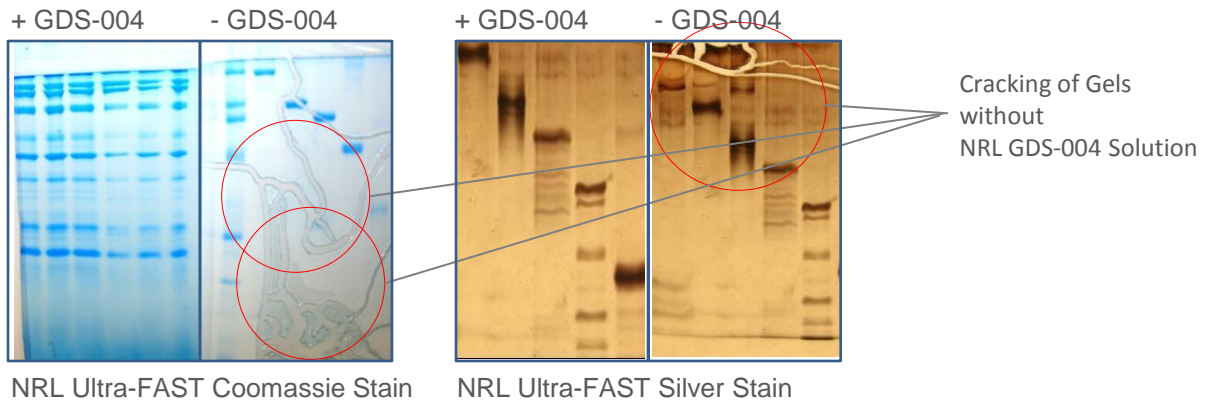

Stability

Dried gels between cellophan membranes are stable for more than 10 years in our laboratories.
 Dried gels can be rehydrated in water at anytime by simply soaking the gel in water for 30 minutes.
 Dried gels can be stored in laboratory note books as original gels.
 Unstained dried gels or gels stained with reversible stain can be rehydrated and probed with antibodies
 Dried gels can be rehydrated and transferred onto nitrocellulose membranes after years of storage.

Ordering Information

Catalog #	Product Description	Price, \$
GDS-001	Easy Dry™ Complete Gel Drying System (1 Unit)	120.00
GDS-002	Easy Dry™ Gel Drying Frame Set (16 x 14 cm)	70.00
GDS-003	Easy Dry™ Gel Drying Frame Sheets(pack of 50) System	70.00
GDS-004	Easy Dry™ Gel Drying Solution (1x 500ml)	75.00
GDS-005	Easy Dry™ Gel Drying Rust-Free Clamps (8 units)	35.00